



OUR EQUIPMENT. YOUR SUCCESS.

TrueDrop™ 500

TECHLAB
SYSTEMS

DROP TEST SYSTEM

Precision drop testing at an affordable price. The TrueDrop™ 500 drop tester. This drop tester represents the ultimate in product and package performance testing and is designed to meet a wide variety of global testing needs.

TEST WITH CONFIDENCE


One of three models in the TRUEDROP™ Series. TRUEDROP™ 500 helps determine the ruggedness of packages and the effectiveness of their interior cushioning when exposed to impact during material handling. The effects of dropping are easily duplicated in a laboratory setting using controlled and precise flat, corner, and edge drop tests.

Complies with ISTA, ASTM, ISO, MIL STD, and other internationally recognized test standards.

Features

Hydraulic Lift System with Multiple-Tine Drop Carriage configured for a variety of products

- 500 lbs (227 kg) maximum payload capacity
- 48" W x 30" D (1219x 762 mm) steel impact surface - Standard size - (custom available)
- 2" - 48" (51- 1200 mm) adjustable drop height Standard drop system - (custom heights available)
- Fully Automated Operator Control Pedestal Station for safe and efficient operation
 - Full featured touch/visual vibrant display
 - Programmable in SAE or Metric units
 - "Arm" and "Drop/Release" indicators for safe operation
 - Programmable drop height, home and drop count input selections
- Versatile package height and width
- Package Holding Fixture holds package in correct orientation for corner and edge impacts

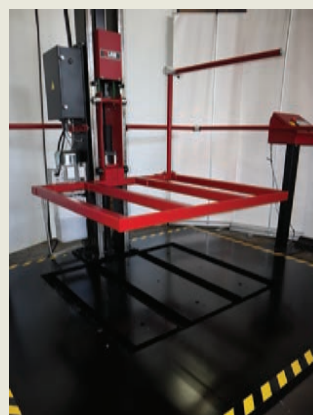

DROP
with CONFIDENCE



Payload Capacity	500 lb (227 kg)/TrueDrop™ 750+ (350 kg)
Drop Height*	2-48 in (51-1200 mm), *higher drop heights available
Maximum Package Size	Dependent on impact surface and drop tines
Impact Zone**	48x30 in (1219x762 mm) - WxD, **expandable as required
Utility Requirement (Electrical)	220-240 VAC/3 Ph/50-60Hz/15 amp or 380-460 VAC/3Ph/50-60Hz/10 amp
Machine Weight	2,000 lb (970 kg)
Machine Dimensions***	96x48x44 in (2.4x2.1x1.1 m) - HxWxD-***Standard model

Due to our continuous commitment to product development, the above specifications and features may be modified without notice.



OUR EQUIPMENT. YOUR SUCCESS.

A Spectral Dynamics Company



L.A.B. Equipment, Inc., 1549 Ardmore Avenue, Itasca, IL 60143, USA

Phone: 630-595-4288 Fax: 630-595-5196 www.labequipment.com

DROP TESTING

Rev B

TrueDrop™ 500 OPTIONS

- Extended drop height
- Enlarged base plate for larger packages
- Extended support tines for wider packages



**Intuitive and Operator Safe
multi-position control station**



Extended Support Tines



**Expanded Base Plate
for larger test items**

